

(19) World Intellectual Property
Organization
International Bureau



(43) International Publication Date
29 September 2005 (29.09.2005)

PCT

(10) International Publication Number
WO 2005/090945 A1

(51) International Patent Classification⁷: **G01N 15/14**,
21/64, 33/50

(21) International Application Number:
PCT/US2005/008870

(22) International Filing Date: 16 March 2005 (16.03.2005)

(25) Filing Language: English

(26) Publication Language: English

(30) Priority Data:
60/553,502 16 March 2004 (16.03.2004) US

(71) Applicant (for all designated States except US): **AMNIS CORPORATION** [US/US]; 2505 Third Avenue, Suite 210, Seattle, WA 98121 (US).

(72) Inventors; and

(75) Inventors/Applicants (for US only): **GEORGE, Thaddeus, C.** [US/US]; 1712 NW 63rd Street, Seattle, WA 98107 (US). **BASILJI, David, A.** [US/US]; 6538 Greenwood Avenue N., Seattle, WA 98103 (US). **HALL, Brian, E.** [US/US]; 1121 17th Avenue, Apt. A, Seattle, WA 98122 (US). **ORTYN, William, E.** [US/US]; 11546 Matsui

Place NE, Bainbridge Island, WA 98110 (US). **SEO, Michael, J.** [US/US]; 8120 SE 77th Place, Mercer Island, WA 98040 (US). **MORRISSEY, Philip, J.** [US/US]; 4412 153rd Avenue SE, Bellevue, WA 98006 (US). **ZIMMERMAN, Cathleen, A.** [US/US]; 7085 NE Bay Hill Road, Bainbridge Island, WA 98110 (US).


(74) Agents: **HERMANN, Karl, R.** et al.; Seed Intellectual Property Law Group PLLC, Suite 6300, 701 Fifth Avenue, Seattle, WA 98104-7092 (US).

(81) Designated States (unless otherwise indicated, for every kind of national protection available): AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BW, BY, BZ, CA, CH, CN, CO, CR, CU, CZ, DE, DK, DM, DZ, EC, EE, EG, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, MZ, NA, NI, NO, NZ, OM, PG, PH, PL, PT, RO, RU, SC, SD, SE, SG, SK, SL, SM, SY, TJ, TM, TN, TR, TT, TZ, UA, UG, US, UZ, VC, VN, YU, ZA, ZM, ZW.

(84) Designated States (unless otherwise indicated, for every kind of regional protection available): ARIPO (BW, GH,

[Continued on next page]

(54) Title: METHOD FOR IMAGING AND DIFFERENTIAL ANALYSIS OF CELLS

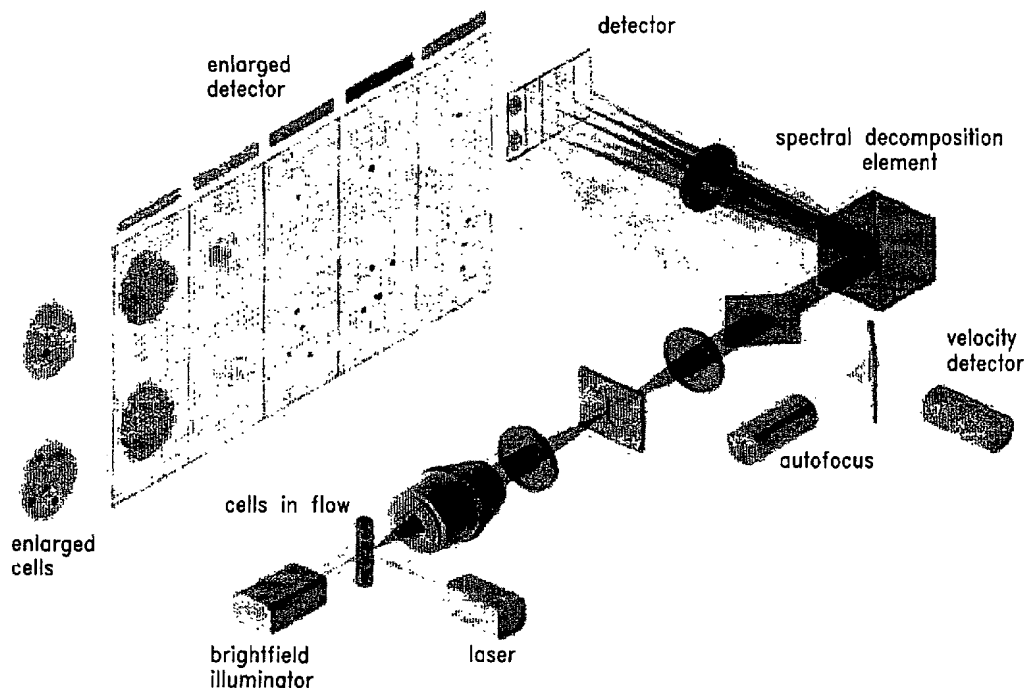

(57) Abstract: Provided are methods for determining and analyzing photometric and morphogenic features of small objects, such as cells to, for example, identify different cell states. In particularly, methods are provided for identifying apoptotic cells, and for distinguishing between cells undergoing apoptosis versus necrosis.

WO 2005/090945 A1



GM, KE, LS, MW, MZ, NA, SD, SL, SZ, TZ, UG, ZM, ZW), Eurasian (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European (AT, BE, BG, CH, CY, CZ, DE, DK, EE, ES, FI, FR, GB, GR, HU, IE, IS, IT, LT, LU, MC, NL, PL, PT, RO, SE, SI, SK, TR), OAPI (BF, BJ, CF, CG, CI, CM, GA, GN, GQ, GW, ML, MR, NE, SN, TD, TG).

Declarations under Rule 4.17:

- as to applicant's entitlement to apply for and be granted a patent (Rule 4.17(ii)) for the following designations AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BW, BY, BZ, CA, CH, CN, CO, CR, CU, CZ, DE, DK, DM, DZ, EC, EE, EG, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, MZ, NA, NI, NO, NZ, OM, PG, PH, PL, PT, RO, RU, SC, SD, SE, SG, SK, SL, SM, SY, TJ, TM, TN, TR, TT, TZ, UA, UG, UZ, VC, VN, YU, ZA, ZM, ZW, ARIPO patent (BW, GH, GM, KE, LS, MW, MZ, NA, SD, SL, SZ, TZ, UG, ZM, ZW), Eurasian patent (AM,

AZ, BY, KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, BG, CH, CY, CZ, DE, DK, EE, ES, FI, FR, GB, GR, HU, IE, IS, IT, LT, LU, MC, NL, PL, PT, RO, SE, SI, SK, TR), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, GQ, GW, ML, MR, NE, SN, TD, TG)

- as to the applicant's entitlement to claim the priority of the earlier application (Rule 4.17(iii)) for all designations
- of inventorship (Rule 4.17(iv)) for US only

Published:

- with international search report
- before the expiration of the time limit for amending the claims and to be republished in the event of receipt of amendments

For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.

METHOD FOR IMAGING AND DIFFERENTIAL ANALYSIS OF CELLS

BACKGROUND OF THE INVENTION

Field of the Invention

The present disclosure relates generally to imaging small moving objects or particles to analyze and distinguish such objects, and more specifically, to a system and method for determining and analyzing photometric and morphogenic features of small objects, such as cells to, for example, identify different cell states.

Description of the Related Art

Apoptosis is a complex, tightly regulated process by which a cell orchestrates its own destruction in response to specific internal or external triggers (Jacobson *et al.*, *Cell* 88:347, 1997; Rathmell and Thompson, *Cell* 109 (Supp):S97, 2002), and proceeds in a manner that is designed to prevent damage to surrounding cells and tissues. Apoptotic cells typically appear shrunken, with condensed chromatin and fragmented nuclei. Although plasma membrane integrity is initially preserved, in later stages the plasma membrane becomes compromised and the cells shed apoptotic bodies consisting of organelles, cytoplasm and/or nuclear fragments. Apoptotic cells are rapidly phagocytosed and eliminated *in vivo*, thus preventing the induction of inflammatory responses, which is a process critical to the maintenance of tissue and immune cell development and homeostasis (Jacobson *et al.*; Rathmell and Thompson; Vaux and Korsmeyer, *Cell* 96:245, 1999). Inappropriately low apoptotic rates can result in cancer or autoimmune disease, while high rates can result in neurodegenerative disease or immunodeficiency (Ashkenazi and Dixit, *Science* 281:1305, 1998; Thompson, *Science* 267:1456, 1995; Fadeel *et al.*, *Leukemia* 14:1514, 2000). In contrast, necrotic cell death is a largely unregulated process in which the cells generally have intact nuclei with limited chromatin condensation. Cells undergoing necrosis do not induce an early phagocytic response. Instead, the cells swell and rupture, and the release of cellular contents can result in significant local tissue damage and inflammation (Jacobson *et al.*).

Research aimed at cell death regulation has produced a number of methods to identify and quantify apoptotic cells, and to distinguish between cells undergoing apoptosis versus necrosis. Among these, flow cytometry has become a commonly used tool in the identification and quantification of apoptosis. Changes in cell size, shape, and granularity associated with apoptosis can be inferred from scattered laser light (Ormerod *et al.*, *J. Immunol. Methods* 153:57, 1992). Early intracellular events, such as the loss of the mitochondrial inner membrane potential or activation and cleavage of caspases, can also be detected using electro-potential sensitive dyes (Castedo *et al.*, *J. Immunol. Methods* 265:39, 2002; Green and Kroemer, *Trends Cell. Biol.* 8:267, 1998; Green and Reed, *Science* 281:1309, 1998; Kroemer and Reed, *Nat. Med.* 6:513, 2000; Lizard *et al.*, *Cytometry* 21:275, 1995) or fluorogenic substrates (Komoriya *et al.*, *J. Exp. Med.* 191:1819, 2000; Smolewski *et al.*, *J. Immunol. Methods* 265:111, 2002; Lecoeur *et al.*, *J. Immunol. Methods* 265:81, 2002). Another early apoptotic event results in exposure of phosphatidylserine on the outer surface of the plasma membrane, which can be detected by fluorochrome-labeled annexin V (van Engeland *et al.*, *Cytometry* 31:1, 1998; Vermes *et al.*, *J. Immunol. Methods* 184:39, 1995; Koopman *et al.*, *Blood* 84:1415, 1994; Verhoven *et al.*, *J. Exp. Med.* 182:1597, 1995). Apoptotic cells eventually lose the ability to exclude cationic nucleotide-binding dyes and nuclear DNA stains with dyes, such as propidium iodide and 7-aminoactinomycin D (7-AAD) (Lecoeur *et al.*, 2002; Gaforio *et al.*, *Cytometry* 49:8, 2002; Ormerod *et al.*, *Cytometry* 14:595, 1993; Schmid *et al.*, *J. Immunol. Methods* 170:145, 1994; Philpott *et al.*, *Blood* 87:2244, 1996). Other techniques that can be used to identify apoptosis include biochemical identification of the activated proteases (*e.g.*, caspases, PARP), release of mitochondrial cytochrome c, quantification of cellular DNA content, and progressive endonucleolytic cleavage of nuclear DNA (Alnemri *et al.*, *Cell* 87:171, 1996; Kohler *et al.*, *J. Immunol. Methods* 265:97, 2002; Gong *et al.*, *Anal. Biochem.* 218:314, 1994; Gorczyca *et al.*, *Leukemia* 7:659, 1993; Gorczyca *et al.*, *Cancer Res.* 53:1945, 1993).

As noted above, conventional flow cytometric methods do not provide direct morphologic evidence of cell death. Indeed, these techniques usually target molecular changes that are associated with apoptosis, but such changes are not always specific to apoptosis and may also be present in cells undergoing necrotic death (Lecoeur *et al.*, 2002; Lecoeur *et al.*, *Cytometry* 44:65, 2001; Kerr *et al.*, *Br. J. Cancer* 26:239, 1972).

For example, necrotic cells, like advanced (late-stage) apoptotic cells, stain with both annexin V and 7-AAD (Lecoeur *et al.*, 2002; Lecoeur *et al.*, 2001). Thus, visualization of the characteristic morphologic changes associated with apoptosis is still considered to be absolutely necessary for its identification (Jacobson *et al.*; Darzynkiewicz *et al. Cytometry* 5 27:1, 1997). Standard microscopic techniques allow visualization of specific molecular and biochemical changes associated with apoptosis and also morphologic changes that distinguish apoptosis from necrosis. However, these standard techniques also require subjective analysis and time-consuming image viewing, which only allows for processing of relatively limited numbers of cells and, therefore, makes it difficult to attain statistically 10 valid comparisons (Tarnok and Gerstner, *Cytometry* 50:133, 2002).

Thus, the need exists for techniques that can provide the statistical power offered by flow cytometry coupled with the objective assessment capabilities associated with microscopic analysis. For example, interest in the dynamic nature of the living cell and efforts to model cell processes (variously termed “cytomics” or “systems biology”) are 15 powerful drivers for new techniques to acquire ever more comprehensive data from cells and cell populations. The present invention meets such needs, and further provides other related advantages.

BRIEF DESCRIPTION OF THE DRAWINGS

Figure 1 shows a schematic representation of the ImageStream 100™ 20 multispectral imaging cytometer.

Figures 2A – 2C show an analysis of cell death using standard flow cytometry and immunofluorescence microscopy. Untreated Jurkat cells (A), Jurkat cells treated with 1 μ M CPT for 18 hrs (B), and Jurkat cells treated with 0.3% hydrogen peroxide for 1 hour (C) were stained with Alexa Fluor 488 conjugated annexin V and 25 7-AAD. Cells were analyzed either by conventional flow cytometry (using a BD FACSort™) or visualized on slides using a Nikon Eclipse E600 fluorescence microscope equipped with bandpass filters appropriate for Alexa Fluor 488 (535/40 nm) and 7-AAD (630/60 nm) fluorescence. The 2-color dot-plots of annexin V vs. 7-AAD, and the brightfield, combined fluorescence and darkfield microscopic images are shown.

Figures 3A – 3E show flow cytometric imaging of untreated, CPT-treated and peroxide-treated Jurkat cells that were stained with Alexa Fluor 488 conjugated annexin V and 7-AAD. Peroxide-treated cells were also separately stained with HLA class I-PE. After staining, equal cell numbers of the three populations of cells were mixed and analyzed by (A) conventional flow cytometry using a FACSort™; and (B) multispectral imaging of cells in flow using an ImageStream 100™ cytometer. The six channel images of cells from representative members of the double positive (DP), single positive (SP), and double negative (DN) populations identified using the ImageStream 100™ are shown in panels C, D and E, respectively.

Figures 4A – 4F show (A) a laser scatter analysis (forward scatter vs. side scatter) of the “DP gate” population of cells from Figure 3A by way of CellQuest™ software on data obtained from a FACSort™ cytometer; (B) a single color histogram of HLA class I-PE on DP cells from Figure 3A as measured using a FACSort™ cytometer; (C) backgating of HLA class I-PE⁺ cells from Figure 4B onto the scatter histogram of Figure 4A, wherein HLA class I-PE⁺ cells (*i.e.*, peroxide-treated, necrotic) are shown in red and HLA class I-PE⁻ cells (*i.e.*, CPT-treated, apoptotic) are shown in blue; (D) a bivariate plot (scatter histogram) of the “Brightfield Area” vs. the “488 nm Scatter Peak Intensity” produced using IDEAS™ software on data obtained using the ImageStream™ cytometer; (E) a single color histogram of HLA class I-PE on DP cells from Figure 3D as measured using the ImageStream™ cytometer; (F) backgating of HLA class I-PE⁺ cells from Figure 4E onto the scatter histogram of Figure 4D, wherein HLA class I-PE⁺ cells (*i.e.*, peroxide-treated, necrotic) are shown in red and HLA class I-PE⁻ cells (*i.e.*, CPT-treated, apoptotic) are shown in yellow.

Figure 5 shows resolution of live, early and late apoptotic, and necrotic cells using morphometric features based on scatter intensity, brightfield area, and nuclear area. Backgating of the four cell populations that had been identified using alternative criteria confirmed their identity as live cells and early apoptotic cells with the DN and SP cells, respectively, shown in Figure 3B (shown in blue and green, center panel), and as necrotic and late apoptotic cells with cells contained in gates R3 and R4, respectively, shown in Figure 4 (shown in yellow and red, center panel).

DETAILED DESCRIPTION OF THE INVENTION

The instant disclosure relates to the use of both photometric and morphometric features derived from multi-mode imagery of objects (*e.g.*, cells) in flow to discriminate cell states or types, and cell features, in heterogeneous populations of cells, including both non-adherent and adherent cell types. A surprising result of the instant disclosure is the ability to discriminate between different cell states, such as differentiating and identifying live cells, necrotic cells, and cells in both the early and late stages of apoptosis, by using unique combinations of features provided in the ImageStream™ Multispectral Imaging Cytometer and the IDEAS™ data analysis software. Discussed in more detail below are single-step methods for basic and complex morphometric classification of objects in flow, which may be combined with comprehensive multispectral imagery and photometric features to allow, for example, the identification of different cell features and/or cell types or states not feasible with standard flow cytometry.

In the present description, any concentration range, percentage range, or integer range is to be understood to include the value of any integer within the recited range and, when appropriate, fractions thereof (such as one tenth and one hundredth of an integer, etc.), unless otherwise indicated. As used herein, the term "about" means $\pm 15\%$. As used herein, the use of an indefinite article, such as "a" or "an", should be understood to refer to the singular and the plural of a noun or noun phrase (*i.e.*, meaning "one or more" of the enumerated elements or components). The use of the alternative (*e.g.*, "or") should be understood to mean either one, both or any combination thereof of the alternatives.

By way of background, methodologies for simultaneous high speed multispectral imaging in brightfield, darkfield, and four channels of fluorescence of cells in flow were recently developed (*see, e.g.*, U.S. Patent Nos. 6,211,955 and 6,249,341). Figure 1 illustrates an exemplary imaging system (*e.g.*, the ImageStream platform). Cells are hydrodynamically focused into a core stream and orthogonally illuminated for both darkfield and fluorescence imaging. The cells are simultaneously trans-illuminated via a spectrally-limited source (*e.g.*, filtered white light or a light emitting diode) for brightfield imaging. Light is collected from the cells with an imaging objective lens and is projected on a charge-coupled detector (CCD). The optical system has a numeric aperture of 0.75 and the CCD pixel size in object space is 0.5 microns square, allowing high resolution

imaging at event rates of approximately 100 cells per second. Each pixel is digitized with 10 bits of intensity resolution, providing a minimum dynamic range of three decades per pixel. In practice, the spread of signals over multiple pixels results in an effective dynamic range that typically exceeds four decades per image. Additionally, the sensitivity of the
5 CCD can be independently controlled for each multispectral image, resulting in a total of approximately six decades of dynamic range across all the images associated with an object.

Prior to projection on the CCD, the light is passed through a spectral decomposition optical system that directs different spectral bands to different lateral
10 positions across the detector (*see, e.g.*, U.S. Patent No. 6,249,341). With this technique, an image is optically decomposed into a set of 6 sub-images, each corresponding to a different color component and spatially isolated from the remaining sub-images. This process allows for identification and quantitation of signals within the cell by physically separating on the detector signals that may originate from overlapping regions of the cell. Spectral
15 decomposition also allows multimode imaging: the simultaneous detection of brightfield, darkfield, and multiple colors of fluorescence. This is exemplified in Figure 1, which depicts a red brightfield illumination source and the associated transmitted light images in the red detector channel adjacent to fluorescent and scattered light images in the other spectral channels. The process of spectral decomposition occurs during the image
20 formation process rather than via digital image processing of a conventional composite image.

The CCD may be operated using a technique called time-delay-integration (TDI), a specialized detector readout mode that preserves sensitivity and image quality even with fast relative movement between the detector and the objects being imaged. As
25 with any CCD, image photons are converted to photocharges in an array of pixels. However, in TDI operation the photocharges are continuously shifted from pixel to pixel down the detector, parallel to the axis of flow. If the photocharge shift rate is synchronized with the velocity of the flowing cell's image, the effect is similar to physically panning a camera: image streaking is avoided despite signal integration times that are orders of
30 magnitude longer than in conventional flow cytometry. For example, an instrument may operate at a continuous data rate of approximately 30 megapixels per second and integrate

signals from each object for 10 milliseconds, allowing the detection of even faint fluorescent probes within cell images that are acquired at high-speed. Careful attention to pump and fluidic system design to achieve highly laminar, non-pulsatile flow eliminates any cell rotation or lateral translation on the time scale of the imaging process (*see, e.g.,*
5 U.S. Patent No. 6, 532,061).

A real-time algorithm analyzes every pixel read from the CCD to detect the presence of object images and calculate a number of basic morphometric and photometric features, which can be used as criteria for data storage. Data files encompassing 10,000-20,000 cells are typically about 100 MB in size and, therefore, can be stored and
10 analyzed using standard personal computers. The TDI readout process operates continuously without any "dead time", which means every cell can be imaged and the coincidental imaging of two or more cells at a time, as depicted in Figure 1, presents no barrier to data acquisition.

Such an imaging system can be employed to determine morphological,
15 photometric, and spectral characteristics of cells and other objects by measuring optical signals, including light scatter, reflection, absorption, fluorescence, phosphorescence, luminescence, etc. As used herein, morphological parameters may be basic (*e.g.,* nuclear shape) or may be complex (*e.g.,* identifying cytoplasm size as the difference between cell size and nuclear size). For example, morphological parameters may include nuclear area,
20 perimeter, texture or spatial frequency content, centroid position, shape (*i.e.,* round, elliptical, barbell-shaped, etc.), volume, and ratios of any of these parameters. Morphological parameters may also include cytoplasm size, texture or spatial frequency content, volume and the like, of cells. As used herein, photometric measurements with the
25 aforementioned imaging system can enable the determination of nuclear optical density, cytoplasm optical density, background optical density, and the ratios of any of these values. An object being imaged can be stimulated into fluorescence or phosphorescence to emit light, or may be luminescent wherein light is produced without stimulation. In each case, the light from the object may be imaged on a TDI detector of the imaging system to determine the presence and amplitude of the emitted light, the number of discrete positions
30 in a cell or other object from which the light signal(s) originate(s), the relative placement of

the signal sources, and the color (wavelength or waveband) of the light emitted at each position in the object.

The present disclosure provides methods of using both photometric and morphometric features derived from multi-mode imagery of objects in flow. Such methods
5 can be employed as a cell analyzer to determine one or more cell states or types, and cell features, in heterogeneous populations of cells when entrained in a fluid flowing through an imaging system. As used herein, cell states or types may include live cells, cells early or late in the process of dying (*e.g.*, apoptotic cells or necrotic cells), cells propagating (*e.g.*, cells in different phases of division), populations and subpopulations of cells (*e.g.*,
10 leucocyte subpopulations in blood), etc., and combinations thereof. However, it should also be understood that these exemplary methods might be used for imaging and distinguishing other moving objects that have identifiable photometric and morphometric features. As used herein, gating refers to a subset of data relating to photometric or morphometric imaging. For example, a gate may be a numerical or graphical boundary of
15 a subset of data that can be used to define the characteristics of particles to be further analyzed. Here, gates have been defined, for example, as a plot boundary that encompasses viable (normal) cells as double negatives (DN gate), or early apoptotic cells as single positives (SP gate), or late apoptotic and necrotic cells as double positives (DP gate). Further, backgating may be a subset of the subset data. For example, a forward scatter
20 versus a side scatter plot in combination with a histogram from an additional marker (*e.g.*, HLA-class I-PE) may be used to backgate a subset (*e.g.*, late apoptotic cells) within the initial subset (*e.g.*, late apoptotic and necrotic cells).

In using an imaging system as described herein, it should be made clear that a separate light source is not required to produce an image of the object (cell), if the object
25 is luminescent (*i.e.*, if the object produces light). However, many of the applications of an imaging system as described herein will require that one or more light sources be used to provide light that is incident on the object being imaged. A person having ordinary skill in the art will know that the location of the light sources substantially affects the interaction of the incident light with the object and the kind of information that can be obtained from the
30 images on a TDI detector.

In addition to imaging an object with the light that is incident on it, a light source can also be used to stimulate emission of light from the object. For example, a cell having been contacted with probe conjugated to a fluorochrome (*e.g.*, such as FITC, PE, APC, Cy5, or Cy5.5) will fluoresce when excited by light, producing a corresponding characteristic emission spectra from any excited fluorochrome probe that can be imaged on a TDI detector. Light sources may alternatively be used for causing the excitation of fluorochrome probes on an object, enabling a TDI detector to image fluorescent spots produced by the probes on the TDI detector at different locations as a result of the spectral dispersion of the light from the object that is provided by prism. The disposition of these fluorescent spots on the TDI detector surface will depend upon their emission spectra and their location in the object.

Each light source may produce light that can either be coherent, non-coherent, broadband or narrowband light, depending upon the application of the imaging system desired. Thus, a tungsten filament light source can be used for applications in which a narrowband light source is not required. For applications such as stimulating the emission of fluorescence from probes, narrowband laser light is preferred, since it also enables a spectrally decomposed, non-distorted image of the object to be produced from light scattered by the object. This scattered light image will be separately resolved from the fluorescent spots produced on a TDI detector, so long as the emission spectra of any of the spots are at different wavelengths than the wavelength of the laser light. The light source can be either of the continuous wave (CW) or pulsed type, preferably a pulsed laser. If a pulsed type illumination source is employed, the extended integration period associated with TDI detection can allow the integration of signal from multiple pulses. Furthermore, it is not necessary for the light to be pulsed in synchronization with the TDI detector.

In the embodiments of the present invention, it is to be understood that relative movement exists between the object being imaged and the imaging system. In most cases, it will be more convenient to move the object than to move the imaging system. However, it is also contemplated that in some cases, the object may remain stationary and the imaging system move relative to it. As a further alternative, both the imaging system and the object may be in motion, which movement may be in different directions and/or at different rates.

In certain aspects, there is provided a method for identifying a specific cell, comprising directing incident light at a cell, using a detector to obtain a side scatter image, and using the spatial frequency content of the side scatter image to identify a specific cell. Within certain embodiments, the methods of the instant disclosure may be used to identify
5 a specific cell subpopulation that is part of larger heterogeneous cell population. For example, the methods of this disclosure may be used to identify a normal cell, a cell undergoing apoptosis (including early and late stage apoptosis), and a cell undergoing necrosis. Alternatively, the methods of the instant disclosure may be used to identify cells at particular stages of replication (S phase, G phase, M phase, etc.). Thus, in a
10 heterogeneous population of cells, the methods of the invention may be used to identify at least one apoptotic cell and at least one necrotic cell and at least one normal (viable) cell. In addition, early stage and late stage apoptotic cells may be identified.

In another aspect, the instant disclosure provides a method for identifying a specific cell, comprising directing incident light at a cell, using a detector to obtain a
15 brightfield image, and using the spatial frequency content of the brightfield image to identify a specific cell. In certain embodiments, the spatial frequency content analyzed is of the nucleus. Any of the aforementioned embodiments may be used within the context of this aspect of the invention.

In a further aspect, the instant disclosure provides a method for identifying a
20 specific cell, comprising contacting a cell with a nuclear marker, directing incident light at the marked cell, using a detector to obtain an image of the cell, and using the nuclear marker image in combination with the spatial frequency content of the cell image to identify a specific cell. Again, any of the previous embodiments may be used within this method. In certain embodiments, only a single nuclear marker is used, such as 7-AAD.

25 In any of the aforementioned methods, multiple images may be collected simultaneously. Furthermore, in any of the aforementioned methods, there is relative motion between the cell and the detector. In addition, in any of the aforementioned methods, the detector is a time delay integration charge-coupled detector.

The instant disclosure also provides a kit for use in a multispectral imaging
30 system to identify a specific cell type, comprising a single nuclear marker, wherein a cell

contacted with the single marker for a time sufficient to allow identification of an apoptotic cell or a necrotic cell with the multispectral imaging system, as described herein.

All of the U.S. patents, U.S. patent application publications, U.S. patent applications, foreign patents, foreign patent applications and non-patent publications
5 referred to in this specification and/or listed in the Application Data Sheet, are incorporated herein by reference, in their entirety. The invention having been described, the following examples are intended to illustrate, and not limit, the invention.

EXAMPLES

EXAMPLE 1

INDUCTION OF APOPTOSIS

Human acute T leukemic Jurkat cell line was obtained from ATCC
5 (Rockville, MD; catalog number CRL-1990) and maintained in RPMI 1640 (Gibco, Grand
Island, NY) containing 5% fetal bovine serum, 1 mM sodium pyruvate (Mediatech,
Herndon, VA), 100 μ M nonessential amino acids, 100 U/ml penicillin, 100 μ g/ml
streptomycin, and 2 mM L-glutamine (BioWhittaker, Walkersville, MD) in 5% CO₂
atmosphere at 37°C. The density of exponentially growing cells was less than 3×10^5 cells
10 per ml at the time of all treatments. To induce apoptosis, cells were treated for 18 hours
with 1 μ M camptothecin (CPT, Sigma), a DNA topoisomerase I inhibitor.

EXAMPLE 2

INDUCTION OF NECROSIS

Human acute T leukemic Jurkat cell line was obtained from ATCC
15 (Rockville, MD; catalog number CRL-1990) and maintained in RPMI 1640 (Gibco, Grand
Island, NY) containing 5% fetal bovine serum, 1 mM sodium pyruvate (Mediatech,
Herndon, VA), 100 μ M nonessential amino acids, 100 U/ml penicillin, 100 μ g/ml
streptomycin, and 2 mM L-glutamine (BioWhittaker, Walkersville, MD) in 5% CO₂
atmosphere at 37°C. The density of exponentially growing cells was less than 3×10^5 cells
20 per ml at the time of all treatments. To induce necrosis, cells were treated for 1 hour with
0.3% hydrogen peroxide (Sigma, St. Louis, MO).

EXAMPLE 3

STAINING TO IDENTIFY APOPTOTIC CELLS AND NECROTIC CELLS

Control (untreated) cell, CPT treated (apoptotic) cells, and peroxide treated
25 (necrotic) cells were independently counted and washed once in phosphate buffered saline
(PBS, Fair Lawn, NJ). Each cell group was resuspended at 10^7 cells/ml in annexin V
Binding Buffer (BD Pharmingen, San Diego, CA) containing Alexa Fluor 488 annexin V
(Molecular Probes, Eugene, Oregon) and 10 μ M 7-aminoactinomycin D (7-AAD,

Molecular Probes) for 10 minutes at room temperature. Necrotic cells were additionally stained with phycoerythrin (PE)-labeled anti-HLA-A, B, C (clone G46-2.6, BD Pharmingen; anti-HLA class I). Each cell group was washed in annexin V Binding Buffer, fixed in 2% paraformaldehyde (Sigma), and analyzed as either single populations or as a mixture by flow cytometry and immunofluorescence microscopy.

EXAMPLE 4

CONVENTIONAL FLOW CYTOMETRY AND IMAGING FLOW CYTOMETRY

For flow cytometry, cell fluorescence data excited by a 488 nm laser were acquired using the FACSort™ cytometer (BD Immunocytometry Systems, San Jose, CA) and analyzed using CellQuest™ (BD Immunocytometry Systems). For imaging flow cytometry, fixed cells at 5×10^7 cells per ml were run at 100 cells per second on an ImageStream100™ (“Beta” version), and the data analyzed using the ImageStream Data Analysis and Exploration Software™ (IDEAS™).

EXAMPLE 5

INSTRUMENTATION FOR MULTISPECTRAL IMAGING FLOW CYTOMETRY

Figure 1 provides an exemplary layout of the ImageStream™ platform. Cells are hydrodynamically focused into a core stream and orthogonally illuminated for both darkfield and fluorescence imaging. The cells are simultaneously trans-illuminated via a spectrally-limited source (*e.g.*, filtered white light or a light emitting diode) for brightfield imaging. Light is collected from the cells with an imaging objective lens and is projected on a charge-coupled detector (CCD). The optical system has a numeric aperture of 0.75 and the CCD pixel size in object space is 0.5 microns square, allowing high resolution imaging at event rates of approximately 100 cells per second. Each pixel is digitized with 10 bits of intensity resolution, providing a minimum dynamic range of three decades per pixel. In practice, the spread of signals over multiple pixels results in an effective dynamic range that typically exceeds four decades per image. Additionally, the sensitivity of the CCD can be independently controlled for each multispectral image,

resulting in a total of approximately six decades of dynamic range across all the images associated with an object.

Prior to projection on the CCD, the light is passed through a spectral decomposition optical system that directs different spectral bands to different lateral positions across the detector (*see, e.g.*, U.S. Patent No. 6,249,341). With this technique, an image is optically decomposed into a set of 6 sub-images, each corresponding to a different color component and spatially isolated from the remaining sub-images. This is exemplified in Figure 1, which depicts a red brightfield illumination source and the associated transmitted light images in the red detector channel adjacent to fluorescent and scattered light images in the other spectral channels. The process of spectral decomposition occurs during the image formation process rather than via digital image processing of a conventional composite image.

The CCD is operated using time-delay-integration (TDI), in which image photons converted to photocharges in an array of pixels are continuously shifted (at a rate synchronized with the velocity of the flowing cell's image) from pixel to pixel down the detector and parallel to the axis of flow to avoid image streaking. For example, the instrument can operate at a continuous data rate of approximately 30 megapixels per second and integrate signal from each object for 10 milliseconds, which allows the detection of even faint fluorescent probes within cell images that are acquired at high speed. Attention to pump and fluidic system design to achieve highly laminar, non-pulsatile flow can eliminate cell rotation or lateral translation on the time scale of the imaging process (*see, e.g.*, U.S. Patent No. 6, 532,061). Every pixel read from the CCD is analyzed by a real-time algorithm that detects the presence of object images and calculates a number of basic morphometric and photometric features, which can be used as criteria for data storage. Data files encompassing 10,000-20,000 cells can be about 100 MB in size, and are stored and analyzed using standard personal computers.

EXAMPLE 6

IMMUNOFLUORESCENCE MICROSCOPY

Fixed control and treated cells were placed on a conventional glass slide (Erie Scientific, Portsmouth, NH), mixed 1:1 with Antifade (Molecular Probes) and

covered with a cover slip. The cells were visualized at 400X using an Eclipse E600 (Nikon, Melville, NY) fluorescence microscope equipped with filters appropriate for Alexa Fluor 488 (535/40 nm emission) and 7-AAD (630/60 nm emission).

EXAMPLE 7

5 CONVENTIONAL ANALYSIS OF CELLS INDUCED TO UNDERGO APOPTOSIS OR NECROSIS

Jurkat T cells were treated with peroxide (to induce necrosis), CPT (to induce apoptosis, which contained cells in both early and late stages of apoptosis), or were untreated (control). The three cell populations were then stained with Alexa Fluor 488 annexin V and 7-AAD and evaluated by brightfield, darkfield, and fluorescence
10 microscopy, and by conventional flow cytometry (Figure 2). The vast majority (>98%) of the control cells were viable at the time of staining, and were annexin V⁻, 7-AAD⁻ (double negative, DN; Figure 2A). CPT-treated (apoptotic) cells had two populations of cells, those that were annexin V⁺ (single positive, SP, or early apoptotic cells), and those that were annexin V⁺, 7-AAD⁺ (double positive, DP, or late apoptotic cells) (Figure 2B).
15 Similar to late apoptotic cells, peroxide-treated (necrotic cells) also stained positively with both annexin V and 7-AAD (Figure 2C). However, the condensed, fragmented nuclei of late apoptotic cells could be easily distinguished from the intact nuclei of necrotic cells by immunofluorescence microscopy. In addition, apoptotic cells exhibited greater darkfield intensity and texture as compared to necrotic cells (see B and C image panels on right,
20 respectively).

EXAMPLE 8

ANALYSIS OF HETEROGENEOUS CELL POPULATION (NORMAL, APOPTOTIC AND NECROTIC)

A mixture of control, apoptotic, and necrotic Jurkat cells (individually
25 prepared as described in Examples 1 and 2) were analyzed in parallel by conventional flow cytometry and on an ImageStream 100™ (Beta system, multispectral imaging flow cytometer). In this experiment, all cells were stained with Alexa Fluor 488-conjugated annexin V and 7-AAD. Necrotic cells were also stained with PE-conjugated anti-HLA class I before mixing the cell populations to aid in distinguishing necrotic cells from late

stage apoptotic cells, and to permit "backgating" when necessary. On the ImageStream™, each cell was simultaneously imaged in darkfield (488 nm laser side scatter), green fluorescence (500-550 nm, annexin V channel), orange fluorescence (550-600 nm, PE channel), red fluorescence (600-650 nm, 7-AAD channel), and brightfield (660-720 nm).

- 5 Cells were grouped into live (DN), early apoptotic (SP), or double positive (DP, which would include late apoptotic and necrotic cells) populations based on the total intensities of annexin V and 7-AAD staining.

Similar bivariate dot-plots of annexin V and 7-AAD staining were obtained in analyses from both the conventional flow cytometer (*see* Figure 3A) and the
10 multispectral imaging flow cytometer (*see* Figure 3B). However, a unique aspect of data collected on the ImageStream™ is that each data point can be "clicked on" to observe the cell imagery associated with each data point. Consequently, each population gate can be used to perform a "virtual cell sort" by displaying individual images of cells that fall within each gate – for example, representative images of cells contained in the DP, SP and DN
15 gates can be "virtually sorted," as shown in Figures 3C, D and E, respectively, with each image row representing a different cell. Early apoptotic cells (cells in the SP gate) appear slightly shrunken, with more complex brightfield and darkfield morphologies, as compared to live cells in the DN gate. The double positive (DP) population contains cells with two distinct morphologies – one containing small, irregularly shaped cells with condensed,
20 fragmented nuclei; and a second containing larger cells with large, unfragmented nuclei that stained uniformly with 7-AAD. The morphology of these two populations of cells is consistent with cells in the late stage of apoptosis and necrosis, respectively. Thus, in the absence of imagery provided by the multispectral imaging flow cytometer, data obtained from a conventional flow cytometer does not permit discrimination of similarly stained
25 cells, such as late apoptotic cells from necrotic cells.

EXAMPLE 9

CONVENTIONAL METHODS TO DISTINGUISH LATE APOPTOTIC AND NECROTIC CELLS

As noted in Example 8, although advanced apoptotic and necrotic cells differ morphologically, they cannot be distinguished based solely on annexin V and
30 7-AAD fluorescence. Plotting the mixed late apoptotic and necrotic DP population on a

forward scatter vs. side scatter (FSC vs. SSC) plot reveals two distinct populations of cells (*see* Figure 4A). Analysis of the DP population of cells obtained on the conventional flow cytometer for staining with PE (which was used to stain only the necrotic subpopulation of cells) permits separation of the necrotic and apoptotic subpopulations of cells (Figure 4B).

- 5 Backgating the PE positive necrotic population in blue reveals that the low SSC population consists of necrotic cells (Figure 4C). However, without the aid of an extra marker (in this case, anti-HLA class I-PE) or imagery as described in Example 8, data obtained from a conventional flow cytometer does not permit discrimination of similarly stained cells, such as late apoptotic cells from necrotic cells.

10

EXAMPLE 10

MULTISPECTRAL IDENTIFICATION OF LATE APOPTOTIC AND NECROTIC CELLS

- Analysis of the DP population with IDEAS™ software for size (Brightfield Area) and Scatter Peak Intensity also revealed two populations of cells (Figure 4D). The nuclei of cells that fell within the high brightfield area, low scatter peak intensity area (R3) were intact, uniformly stained with 7-AAD, and had a morphology consistent with necrotic cells. The nuclei of cells that fell within the low brightfield area, high scatter peak intensity area (R4) were condensed and fragmented, and had a morphology consistent with cells in the late stages of apoptosis. Backgating PE-positive cells (identified in the histogram shown in Figure 4E) in yellow verified that R3 gated cells were derived from the necrotic treatment group (Figure 4F). This conclusion is further supported by morphologic examination of cells in the image galleries of the R3 and R4 gated cells, and confirms that the low area/high texture cells were apoptotic (HLA-class I PE⁻ cells containing fragmented 7-AAD staining nuclei; lower right gallery) while high area/low texture cells were necrotic (HLA class I-PE⁺ cells containing uniform 7-AAD staining nuclei; upper right gallery).
- 25 Thus, the data obtained from multispectral imaging is provided in form that allows one to distinguish similarly stained cells, such as late apoptotic cells from necrotic cells.

EXAMPLE 11

COMPLEX MORPHOLOGIC FEATURE IDENTITY USING MULTISPECTRAL IMAGING

Multispectral image data collection not only enables calculation of standard intensity-based parameters and statistics employed in conventional flow cytometry, but
5 also permits quantitation of numerous other morphologic features (*e.g.*, cell area, perimeter, aspect ratio, texture, spot counts, cell centroid, gradient intensity, spatial frequency). Using this capability, it is possible to distinguish all four cell populations (*i.e.*, live, early apoptotic, late apoptotic and necrotic cells) in a single step using morphologic features derived from 7-AAD, brightfield and darkfield imagery (and in the absence of other
10 staining procedures often used) to “identify” apoptotic cells.

By subtracting the 7-AAD image area (nuclear size) from the Brightfield area (cell size), a value is obtained that is an indication of cytoplasmic size. When this complex morphologic feature (herein referred to as “Brightfield – 7-AAD Area”) was used in conjunction with a feature derived from darkfield imagery (*i.e.*, 488 nm spatial scatter
15 frequency, which is an indication of internal cell complexity or cell granularity), four subpopulations of cells became evident (*see* Figure 5). The 488 nm spatial scatter frequency can be calculated by computing the standard deviation of the individual pixel intensities within the segmented dark field image mask.

Live cells (depicted in blue, Figure 5, center panel) excluded the
20 cell-impermeant 7-AAD fluorescent DNA binding dye, which minimized the nuclear image area and resulted in cells with a large calculated cytoplasmic area.

Early apoptotic cells (depicted in green, Figure 5, center panel) are just as effective as live cells at excluding 7-AAD, but their total brightfield area is slightly smaller due to the early stages of cytoplasmic blebbing, thereby resulting in an intermediate value
25 for the “Brightfield - 7AAD Area” parameter. Also associated with the early stages of apoptosis is a significant increase in 488 nm scatter peak intensity, which clearly separates these cells from live cells on the vertical axis of the dot-plot.

Necrotic cells (depicted in yellow, Figure 5, center panel) and late apoptotic cells (depicted in red, Figure 5, center panel) both had compromised membrane integrity,
30 which permits free entry of 7-AAD and thus strong nuclear images of relatively large area, shifting these populations to the left on the dot-plot. However, these two cell populations

can be clearly separated based on the peak intensity measurements derived from their 488 nm scatter parameters. Necrotic cells produce darkfield images of relatively low complexity compared to the more complex and heterogeneous darkfield images of late apoptotic cells, thus clearly separating the two populations in the vertical axis.

5 Inspection of the associated image galleries associated with these four gated populations of cells confirmed the classification of each population (*see* Figure 5, upper and lower panels).

10 From the foregoing it will be appreciated that, although specific embodiments of the invention have been described herein for purposes of illustration, various modifications may be made without deviating from the spirit and scope of the invention. Accordingly, the invention is not limited except as by the appended claims.

CLAIMS

1. A method for identifying a specific cell, comprising directing incident light at a cell, using a detector to obtain a side scatter image, and using the spatial frequency content of the side scatter image to identify a specific cell.
2. The method of claim 1 wherein there is relative motion between the cell and the detector.
3. The method of claim 1 wherein a specific cell subpopulation is identified with a heterogeneous cell population.
4. The method of claim 1 wherein the specific cell identified is an apoptotic cell.
5. The method of claim 4 wherein the apoptotic cell is an early stage apoptotic cell or a late stage apoptotic cell.
6. The method of claim 1 wherein the specific cell identified is a necrotic cell.
7. The method of claim 1 wherein the specific cell identified is at least one of an apoptotic cell and a necrotic cell.
8. A method for identifying a specific cell, comprising directing incident light at a cell, using a detector to obtain a brightfield image, and using the spatial frequency content of the brightfield image to identify a specific cell.
9. The method of claim 8 wherein there is relative motion between the cell and the detector.

10. The method of claim 8 wherein a specific cell subpopulation is identified with a heterogeneous cell population.
11. The method of claim 8 wherein the specific cell identified is an apoptotic cell.
12. The method of claim 11 wherein the apoptotic cell is an early stage apoptotic cell or a late stage apoptotic cell.
13. The method of claim 8 wherein the specific cell identified is a necrotic cell.
14. The method of claim 8 wherein the specific cell identified is at least one of an apoptotic cell and a necrotic cell.
15. The method of claim 8 wherein the spatial frequency content is of the nucleus.
16. A method for identifying a specific cell, comprising contacting a cell with a nuclear marker, directing incident light at the marked cell, using a detector to obtain an image of the cell, and using the nuclear marker image in combination with the spatial frequency content of the cell image to identify a specific cell.
17. The method of claim 16 wherein there is relative motion between the cell and the detector.
18. The method of claim 16 wherein a specific cell subpopulation is identified with a heterogeneous cell population.
19. The method of claim 16 wherein the specific cell identified is an apoptotic cell.

20. The method of claim 19 wherein the apoptotic cell is an early stage apoptotic cell or a late stage apoptotic cell.
21. The method of claim 16 wherein the specific cell identified is a necrotic cell.
22. The method of claim 16 wherein the specific cell identified is at least one of an apoptotic cell and a necrotic cell.
23. The method of claim 16 wherein a single nuclear marker is used.
24. The method of claim 16 wherein the single nuclear marker is 7-aminoactinomycin D.
25. The method according to any one of claims 16-24 wherein the images are collected simultaneously.
26. The method according to any one of claims 1-24 wherein the detector is a time delay integration charge-coupled detector.
27. A kit for use in a multispectral imaging system to identify a specific cell type, comprising a single nuclear marker, wherein a cell contacted with the single marker for a time sufficient to allow identification of an apoptotic cell or a necrotic cell with the multispectral imaging system.
28. The kit of claim 27 wherein the single nuclear marker is 7-aminoactinomycin D.

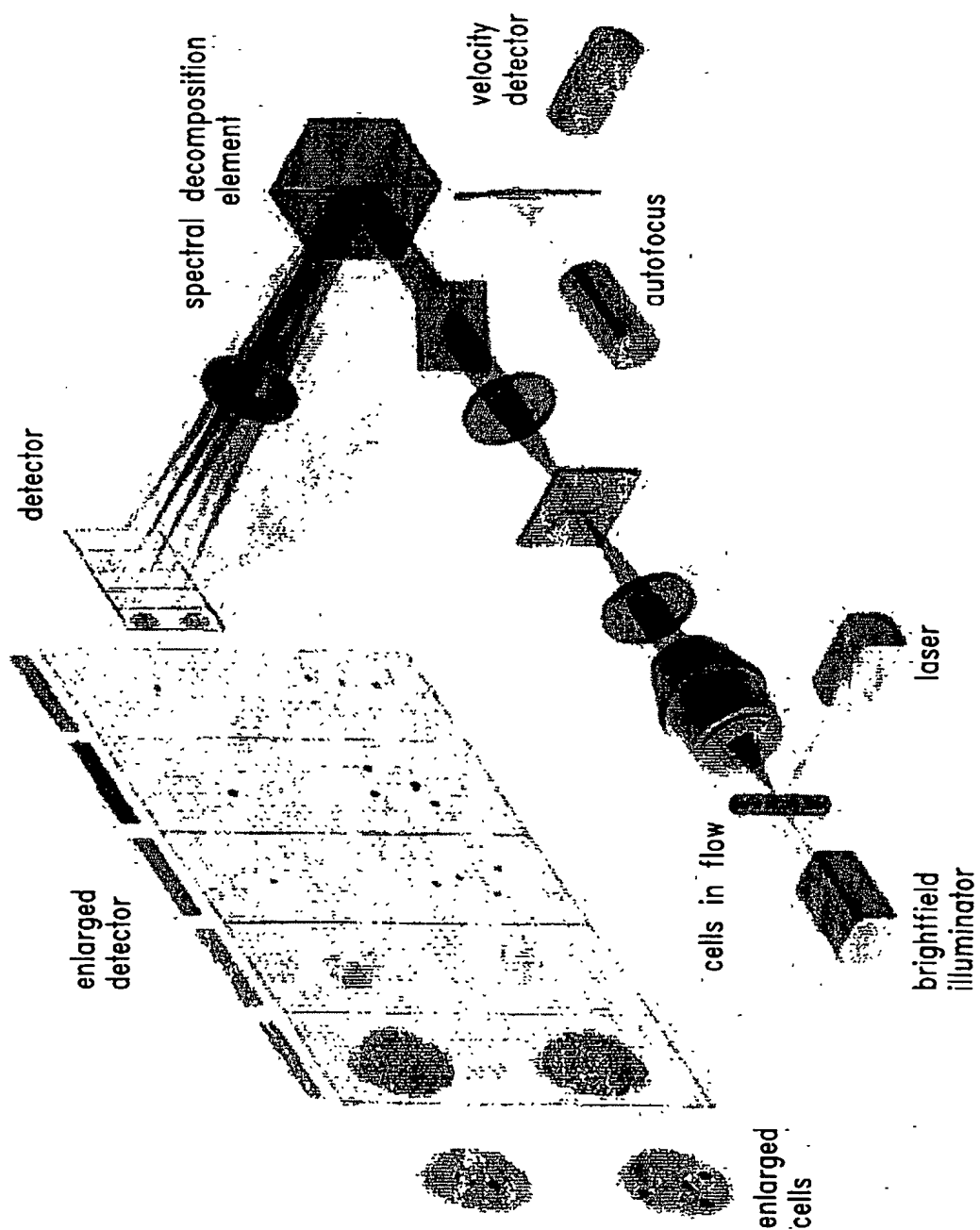


FIG. 1

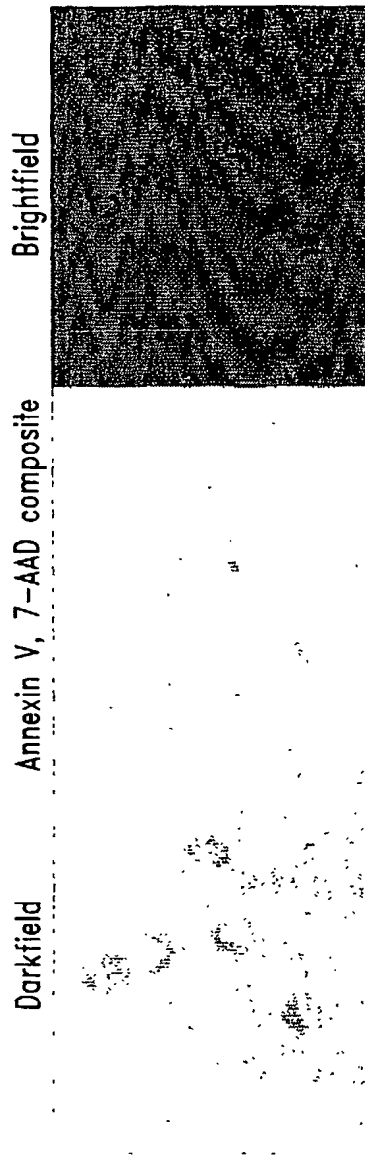


FIG. 2A



FIG. 2B

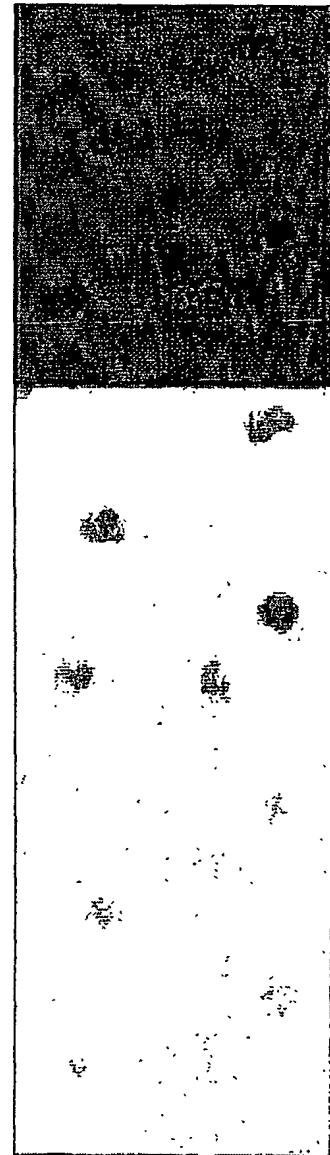
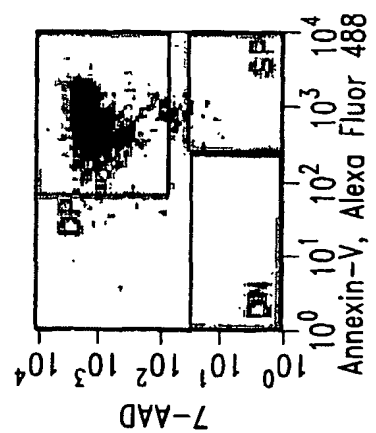
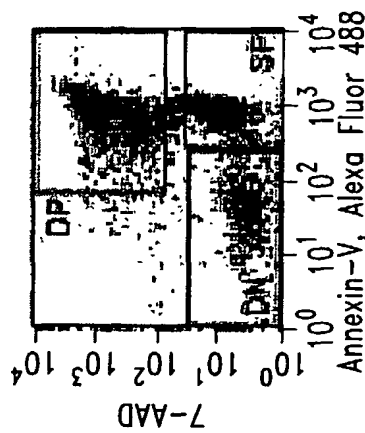
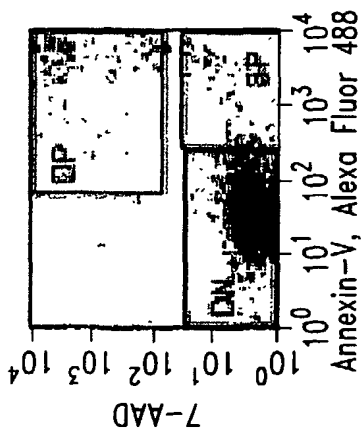


FIG. 2C



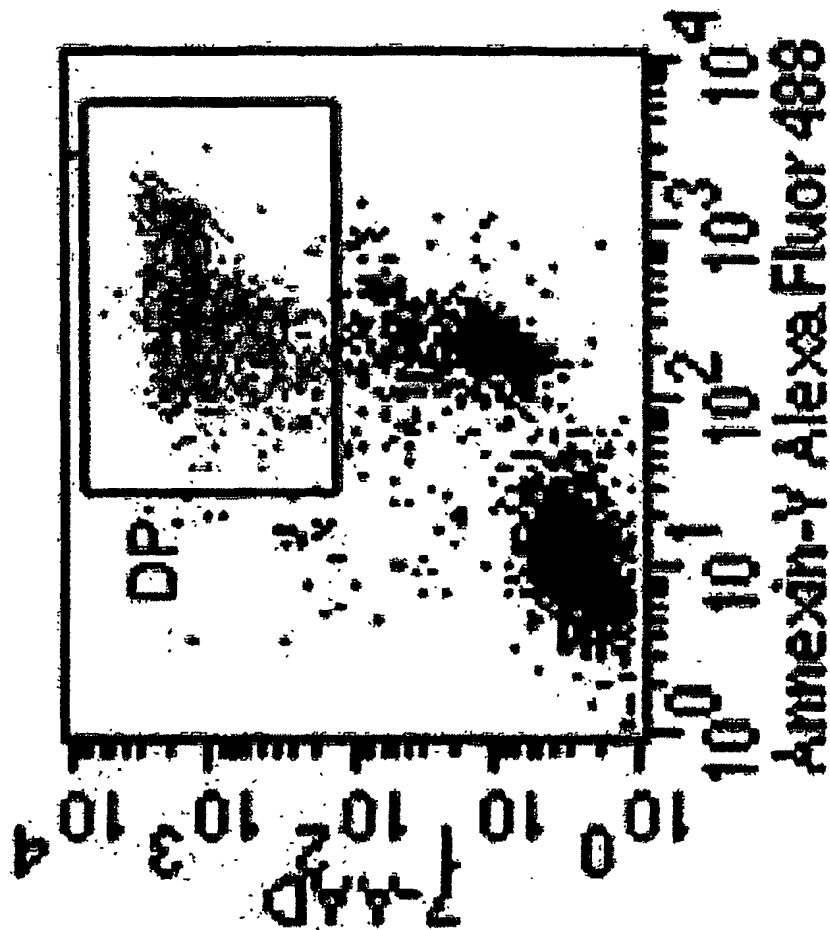


FIG. 3A

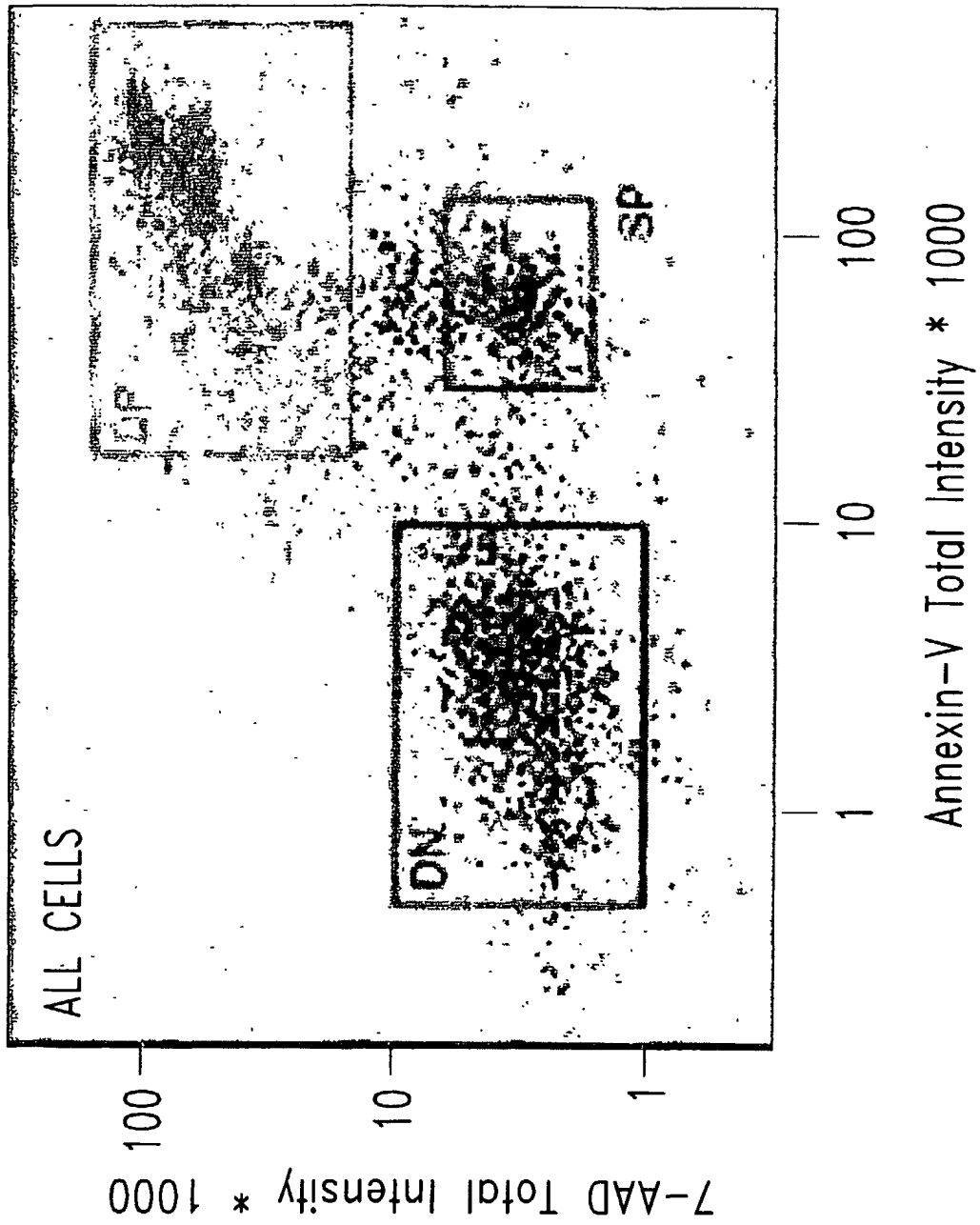


FIG. 3B

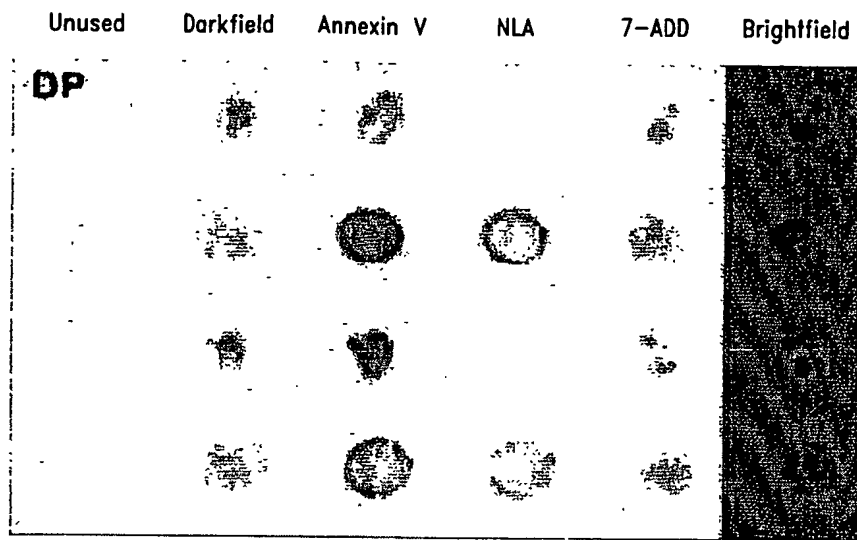


FIG. 3C

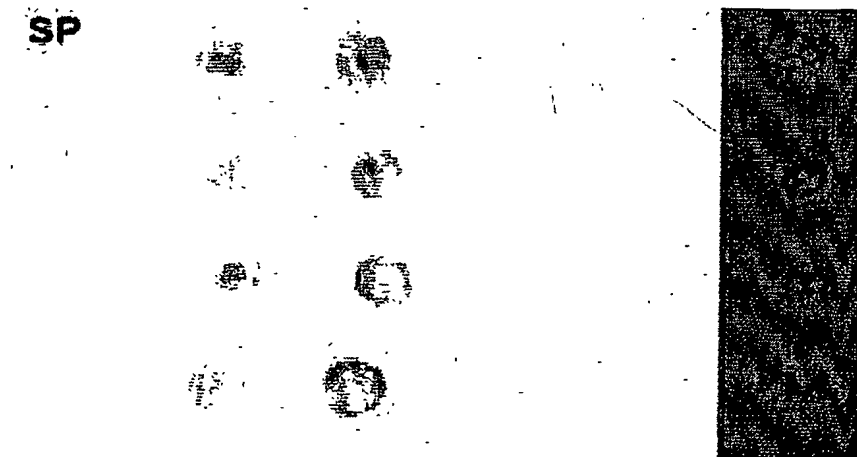


FIG. 3D

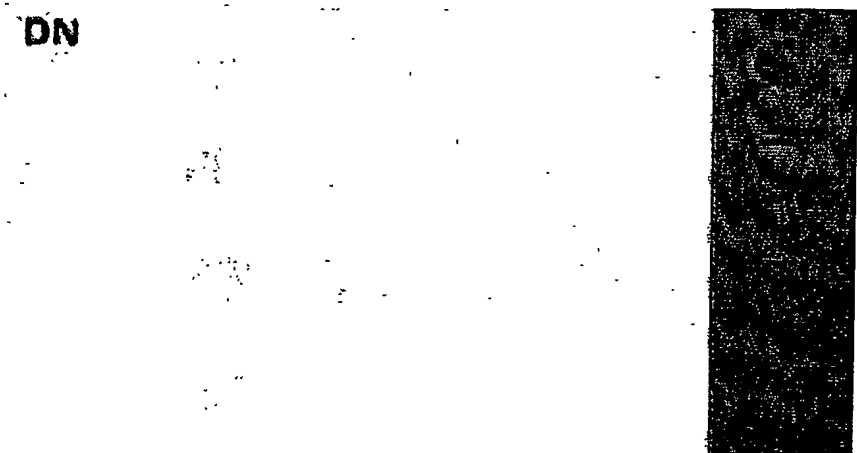
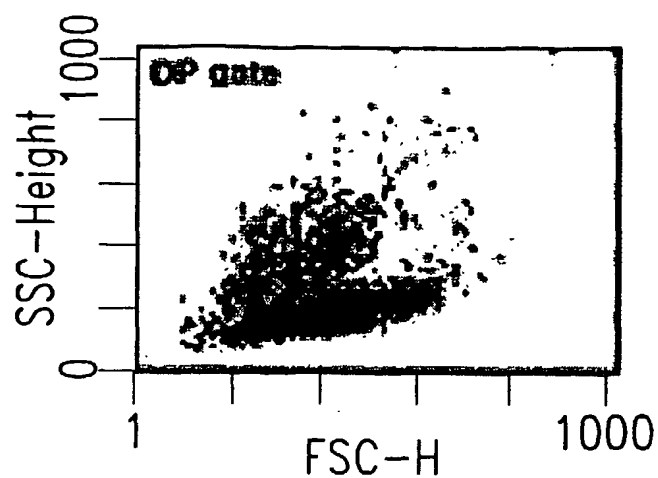
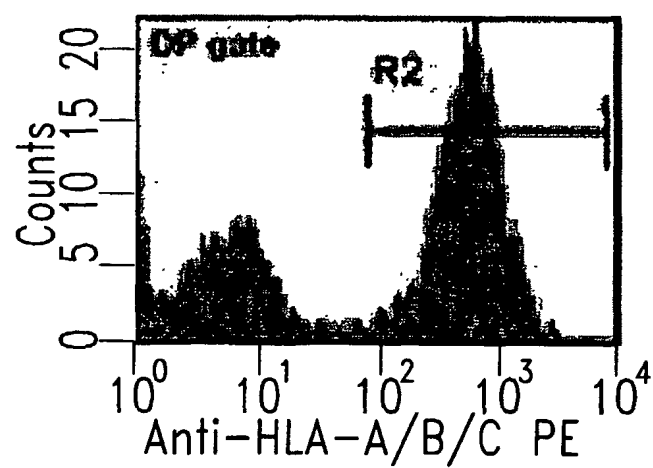
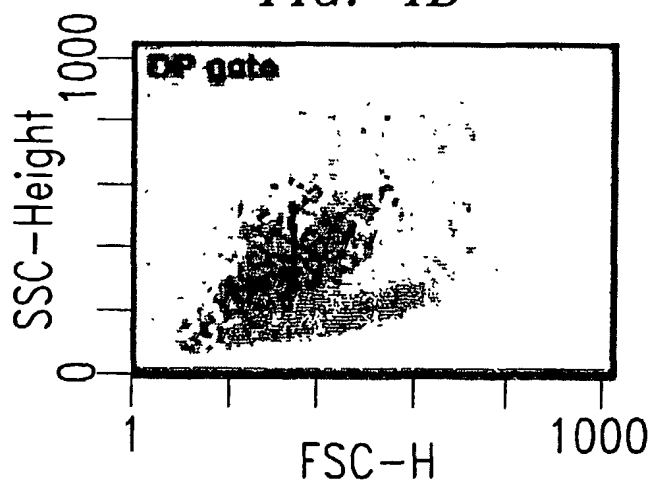


FIG. 3E

*FIG. 4A**FIG. 4B**FIG. 4C*

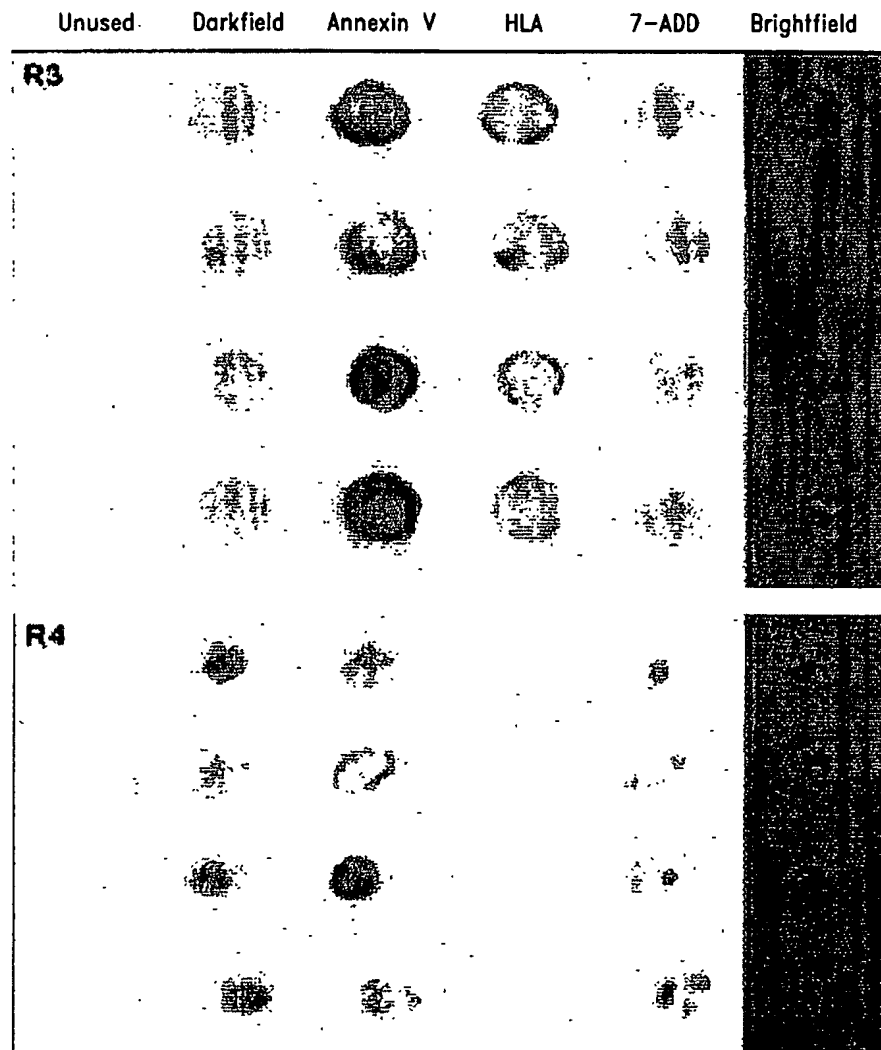
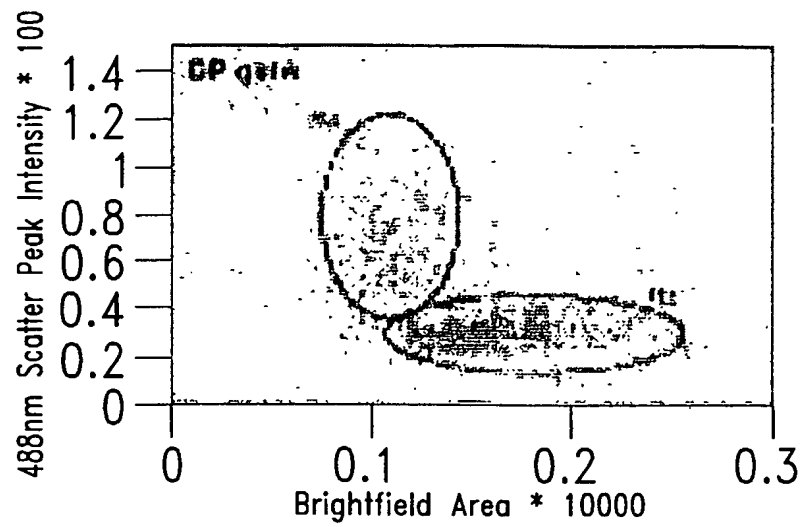
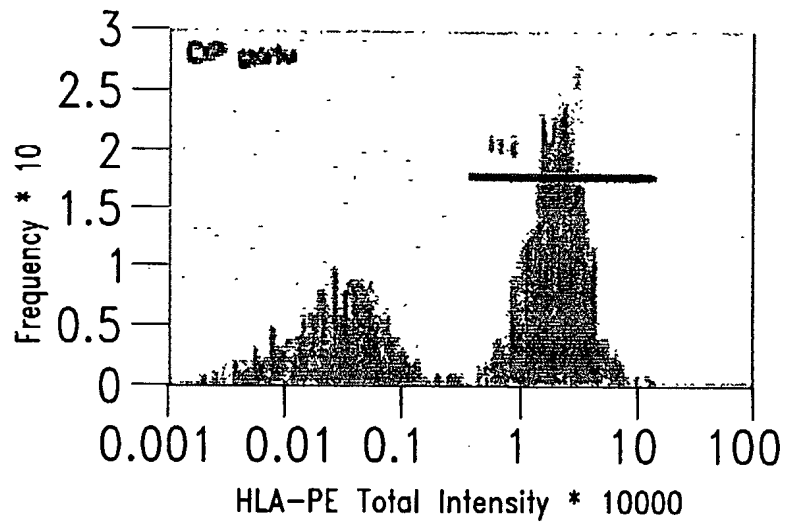
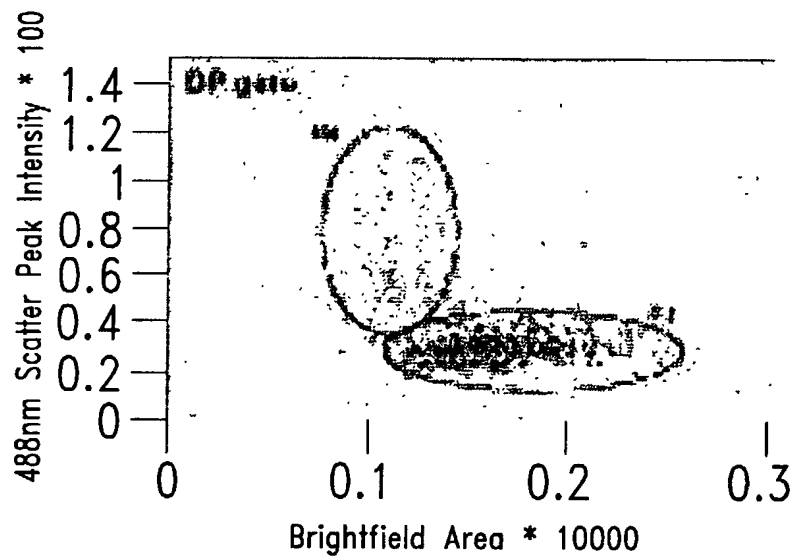
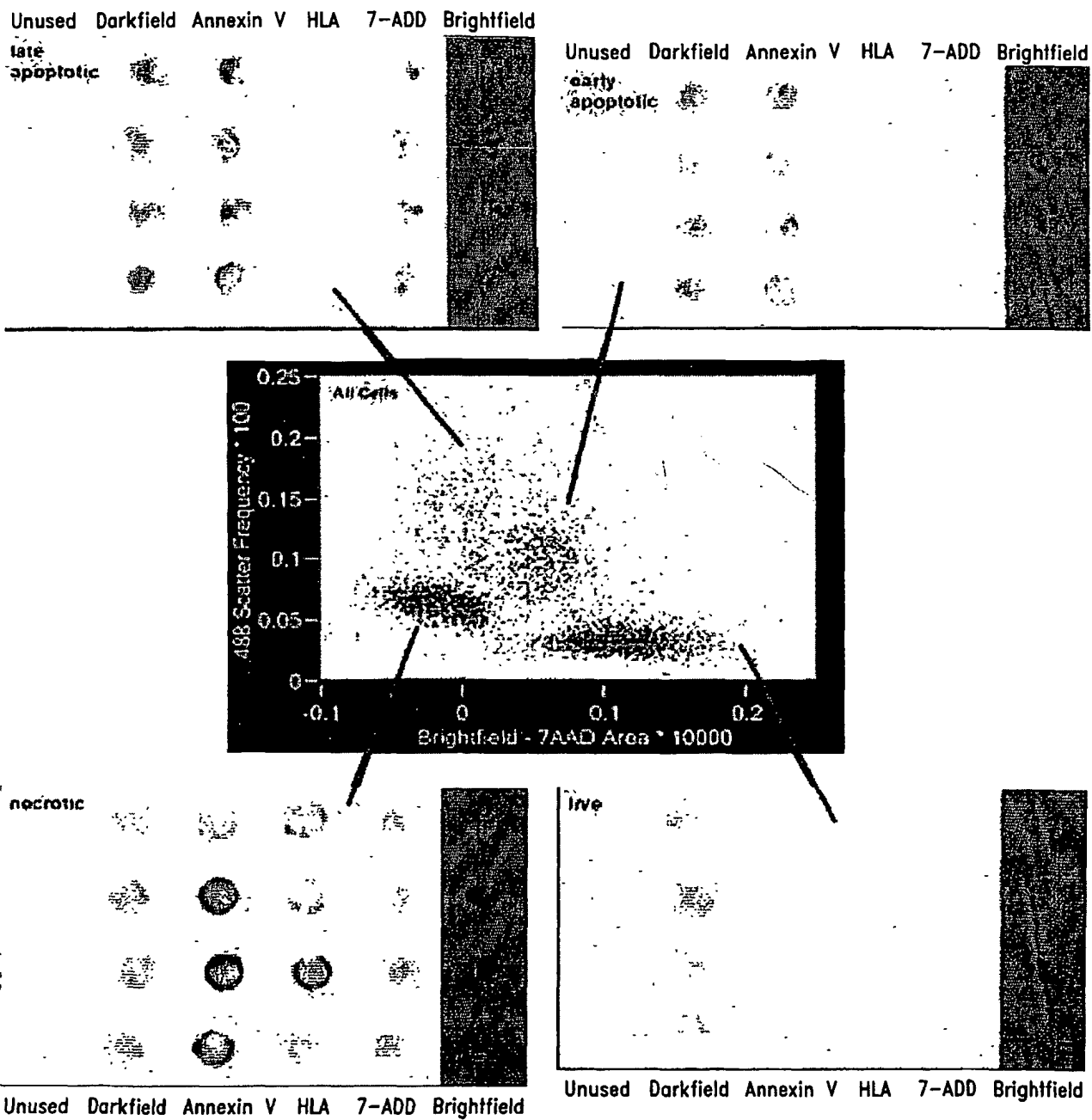


FIG. 4D

*FIG. 4E**FIG. 4F*

**FIG. 5**

INTERNATIONAL SEARCH REPORT

International Application No
PCT/US2005/008870

A. CLASSIFICATION OF SUBJECT MATTER

IPC 7 G01N15/14 G01N21/64 G01N33/50

According to International Patent Classification (IPC) or to both national classification and IPC

B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)

IPC 7 G01N

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practical, search terms used)

EPO-Internal, WPI Data, INSPEC

C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
P, X	GEORGE THADDEUS C ET AL: "Distinguishing modes of cell death using the ImageStream multispectral imaging flow cytometer." CYTOMETRY. PART A : THE JOURNAL OF THE INTERNATIONAL SOCIETY FOR ANALYTICAL CYTOLOGY. JUN 2004, vol. 59, no. 2, June 2004 (2004-06), pages 237-245, XP002334770 ISSN: 1552-4922 the whole document	1-28
X	US 2004/021868 A1 (ORTYN WILLIAM E ET AL) 5 February 2004 (2004-02-05)	1-23, 25-27
Y	paragraphs '0011!, '0063! - '0069!, '0078! - '0080!, '0107!, '0121!, '0135! - '0138!; figure 5	24, 28
	----- -/--	

☒ Further documents are listed in the continuation of box C.

☒ Patent family members are listed in annex.

* Special categories of cited documents :

- *A* document defining the general state of the art which is not considered to be of particular relevance
- *E* earlier document but published on or after the international filing date
- *L* document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)
- *O* document referring to an oral disclosure, use, exhibition or other means
- *P* document published prior to the international filing date but later than the priority date claimed

- *T* later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention
- *X* document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone
- *Y* document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.
- *&* document member of the same patent family

Date of the actual completion of the international search

5 July 2005

Date of mailing of the international search report

19/07/2005

Name and mailing address of the ISA

European Patent Office, P.B. 5818 Patentlaan 2
NL - 2280 HV Rijswijk
Tel. (+31-70) 340-2040, Tx. 31 651 epo nl,
Fax: (+31-70) 340-3016

Authorized officer

Duijs, E

INTERNATIONAL SEARCH REPORT

International Application No
PCT/US2005/008870

C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT		
Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	EP 1 316 793 A (LEIST, CHRISTIAN; BALZER, ANDREAS) 4 June 2003 (2003-06-04)	8,10-16, 18-23,27
Y	paragraphs '0005! - '0011!, '0029! - '0032!, '0050!, '0055!, '0062! - '0064!; figures 1,4 -----	24,28
Y	WO 95/20148 A (COULTER CORPORATION) 27 July 1995 (1995-07-27)	24,28
A	page 7, line 31 - line 37 page 2, line 29 - line 33 -----	3,10,18
A	MAJNO G ET AL: "APOPTOSIS ONCOSIS AND NECROSIS AN OVERVIEW OF CELL DEATH" AMERICAN JOURNAL OF PATHOLOGY, PHILADELPHIA, PA, US, vol. 146, no. 1, January 1995 (1995-01), pages 3-15, XP000884170 ISSN: 0002-9440 the whole document -----	

INTERNATIONAL SEARCH REPORT

Information on patent family members

International Application No

PCT/US2005/008870

Patent document cited in search report		Publication date	Patent family member(s)	Publication date
US 2004021868	A1	05-02-2004	US 2002071121 A1	13-06-2002
			US 2001021018 A1	13-09-2001
			US 6211955 B1	03-04-2001
			US 6249341 B1	19-06-2001
			US 2004218184 A1	04-11-2004
			EP 1334338 A2	13-08-2003
			JP 2004522163 T	22-07-2004
			WO 02101339 A2	19-12-2002
			US 2002030812 A1	14-03-2002
			AU 1199401 A	31-07-2001
			CA 2401614 A1	26-07-2001
			EP 1272822 A1	08-01-2003
			WO 0153783 A1	26-07-2001
			AU 772331 B2	22-04-2004
			AU 1196901 A	31-07-2001
			CA 2395627 A1	26-07-2001
			EP 1257794 A1	20-11-2002
			JP 2003520954 T	08-07-2003
			WO 0153784 A1	26-07-2001
			US 2003137661 A1	24-07-2003
			US 6580504 B1	17-06-2003
			US 2004080748 A1	29-04-2004
			AU 1189802 A	22-04-2002
			AU 1189902 A	22-04-2002
			AU 1191302 A	22-04-2002
			AU 1315702 A	22-04-2002
			EP 1330650 A2	30-07-2003
			JP 2004532969 T	28-10-2004
			WO 0231583 A1	18-04-2002
			WO 0231182 A2	18-04-2002
			WO 0231501 A1	18-04-2002
			WO 0231467 A1	18-04-2002
			US 2002094116 A1	18-07-2002
			US 2002127603 A1	12-09-2002
			US 2002057432 A1	16-05-2002
			US 2002051070 A1	02-05-2002
EP 1316793	A	04-06-2003	EP 1316793 A1	04-06-2003
WO 9520148	A	27-07-1995	AU 1566995 A	08-08-1995
			EP 0745217 A1	04-12-1996
			JP 9508270 T	26-08-1997
			WO 9520148 A1	27-07-1995